


# NETCELL

High Quality ENT PVA Sponge Products



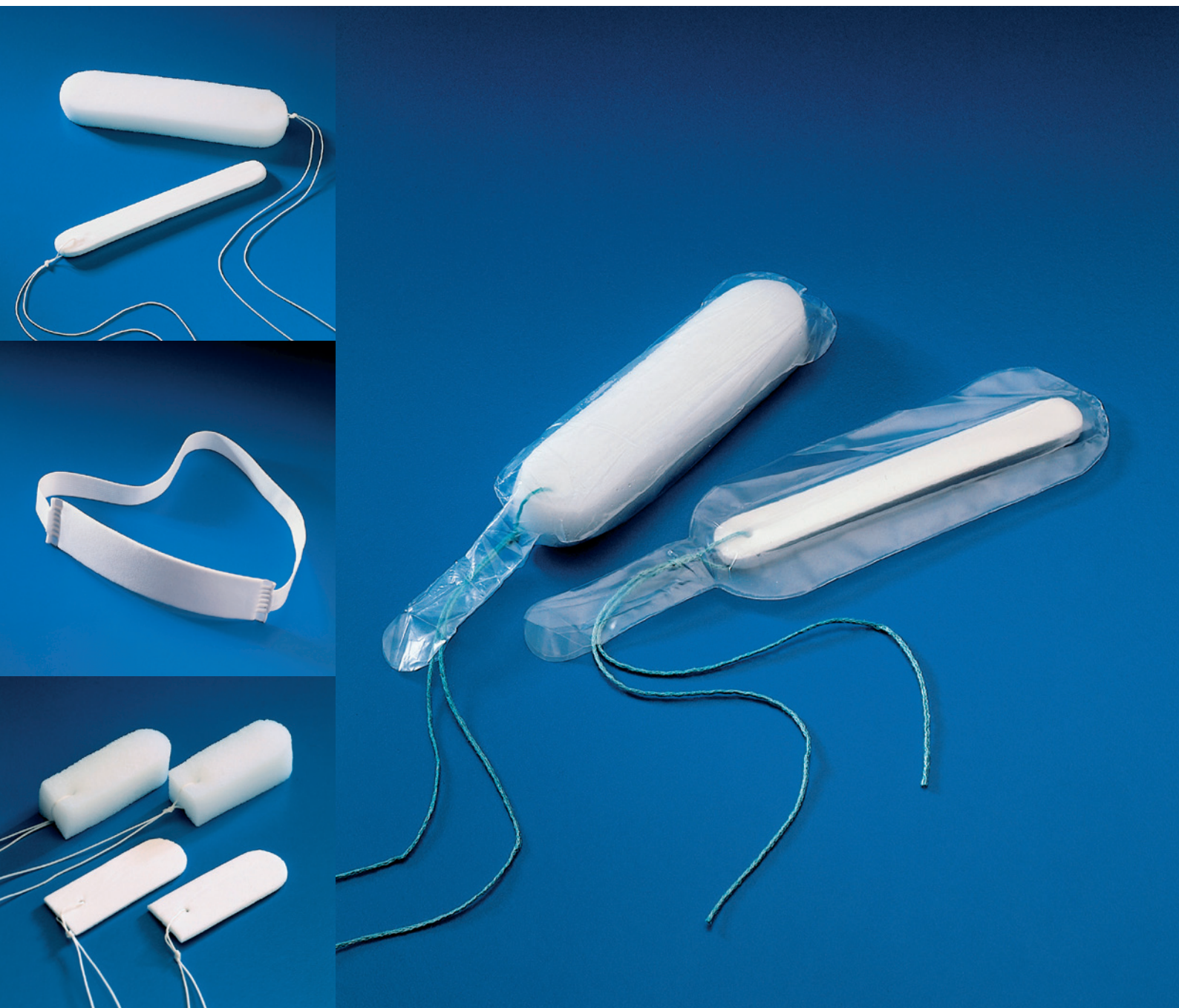

Innovia Medical is a **family of specialty surgical companies** that are passionately committed to partnering with medical professionals to help elevate the delivery of patient care and improve clinical outcomes.

## ELEVATING CARE. IMPROVING OUTCOMES.

---

Innovia Medical companies develop and manufacture a portfolio of product solutions designed with both the medical professional and their patients in mind, to deliver care that:

- Saves time for both the medical professional and patient
- Delivers a consistent and predictable care experience
- Improves patient confidence with reduced recovery times and less discomfort

Innovia Medical does this through a commitment to innovation leadership, unsurpassed product quality and a dedication to provide responsive, personalized customer service.

## INNOVATION LEADERSHIP.

---

Innovation leadership is fueled through our dedication to collaboration with medical professionals. We regularly seek feedback and input from physicians, nurses, and other medical professionals through product research, testing and development.

This allows us to truly understand the solutions medical professionals need, and guides our innovative culture focused on inventing new, industry-leading solutions and reinventing current solutions that make ordinary products extraordinary, across the specialty surgery spectrum.

## UNSURPASSED PRODUCT QUALITY.

---

Our products deliver proven solutions and dependability which is the result of rigorous quality control and a relentless commitment to ongoing manufacturing improvements. We use only the highest quality materials and subject all products to rigorous inspection protocols. This is why we're proud to have so many customers who have continued to rely on our solutions for more than forty years.

## RESPONSIVE, PERSONALIZED CUSTOMER SERVICE.

---

From our sales representatives to our customer service team, every associate at Innovia Medical strives to deliver an unrivaled customer experience. This begins by ensuring that each customer receives the product that is right for their specific needs. We provide assistance through technical support and live customer service representatives with extensive product knowledge, available throughout the day. We deliver quick and thorough responses, every time—from order and shipping confirmations to delivery date communications and more ensuring continuous satisfaction.

We know you're committed to providing the best patient care experience possible and that's why we commit to being your trusted partner helping you every step of the way. At Innovia Medical, we are not just products, we are part of the patient care solution.

Learn more at [www.innoviamedical.com](http://www.innoviamedical.com).

# Reasons to choose NETCELL PVA

High quality PVA sponge is an easy and effective way to improve patient comfort and reduce trauma.

## NETCELL PVA Sponge


- Smooth and soft surface
- Small uniformed pores
- **Decreased** risk of tissue adhesion and postoperative bleeding

Exclusive NETCELL PVA technology features in all our ENT sponge products to increase absorbency and fluid retention whilst decreasing risk of tissue adhesion.




## Competitor PVA Sponge

- Rough and uneven surface
- Large irregular pores
- **Increased** risk of tissue adhesion and postoperative bleeding

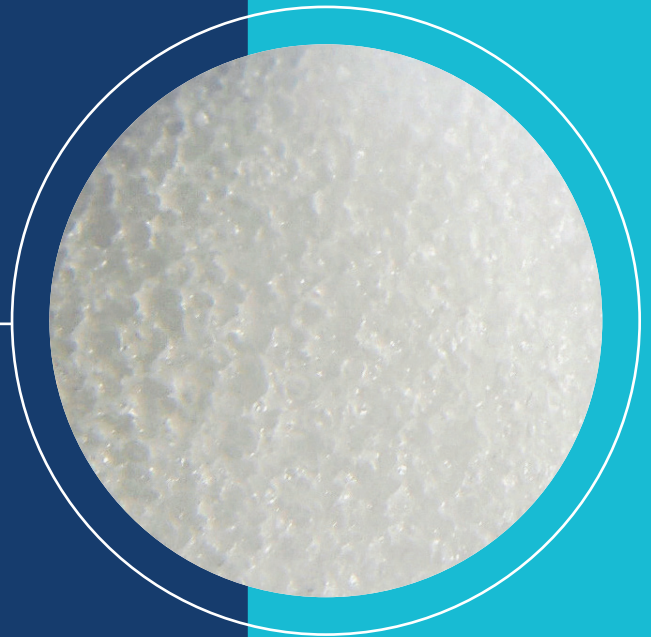



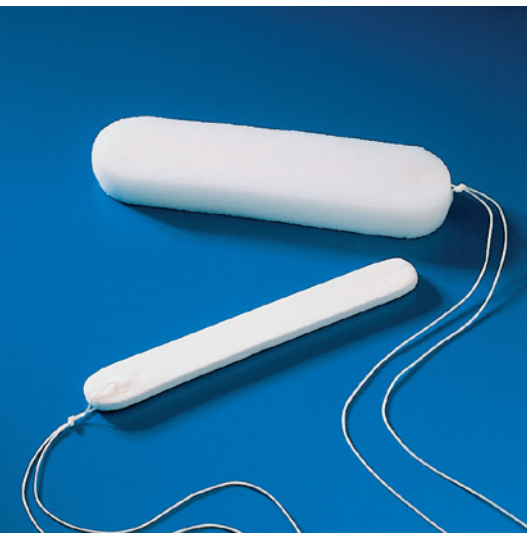

# PVA Nasal Pack

- Postoperative stenting tamponade
- Specifically designed for use after septal, turbinate and rhinoplasty surgery
- Provides gentle, equal pressure within the nasal cavity when absorbing postoperative bleeding
- Available with airway tube for improved venting, patient comfort and reduced risk of obstruction
- Suitable for delivering medication directly to the surgical area

Supplied sterile, single-use

Length: 4.5cm

Description	Dimensions	Pack Size	Code
4.5 cm Nasal Pack	4.5 x 1.5 x 2.0 cm	Box of 20	10-100
4.5 cm Nasal Pack with Tube	4.5 x 1.5 x 2.0 cm	Box of 20	10-100T
4.5 cm Nasal Slimpack	4.5 x 1.0 x 2.0 cm	Box of 20	10-101

Length: 6.0cm

Description	Dimensions	Pack Size	Code
6.0 cm Nasal Pack	6.0 x 1.5 x 2.0 cm	Box of 20	10-105
6.0 cm Nasal Pack with Tube	6.0 x 1.5 x 2.0 cm	Box of 20	10-105T
6.0 cm Nasal Slimpack	6.0 x 1.0 x 2.0 cm	Box of 20	10-106

Length: 8.0cm


Description	Dimensions	Pack Size	Code
8.0 cm Nasal Pack	8.0 x 1.5 x 2.0 cm	Box of 20	10-110
8.0 cm Nasal Pack with Tube	8.0 x 1.5 x 2.0 cm	Box of 20	10-110T
8.0 cm Nasal Slimpack	8.0 x 1.0 x 2.0 cm	Box of 20	10-111
8.0 cm Nasal Slimpack Layered	8.0 x 1.0 x 2.0 cm	Box of 20	10-111L
8.0 cm Nasal Slimpack with Tube	8.0 x 1.0 x 2.0 cm	Box of 20	10-111T

Length: 10.0cm

Description	Dimensions	Pack Size	Code
10.0cm Nasal Pack	10.0 x 1.5 x 2.0 cm	Box of 20	10-115
10.0 cm Nasal Slimpack	10.0 x 1.0 x 2.0 cm	Box of 20	10-116

# PVA Anatomical Nasal Pack

- Postoperative stenting tamponade
- Ergonomically designed to fill the nasal cavity following septal, turbinate or rhinoplasty surgery
- Provides gentle, equal pressure within the nasal cavity when absorbing postoperative bleeding
- Available with airway tube for improved venting, patient comfort and reduced risk of obstruction
- Suitable for delivering medication directly to the surgical area

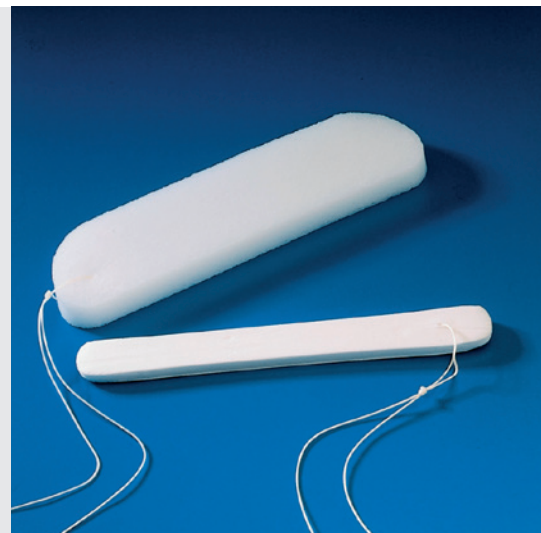


Supplied sterile, single-use

Length: 8.0cm



Description	Dimensions	Pack Size	Code
Anatomical Pack	8.0 x 1.5 x 3.0 cm	Box of 10	10-120
Anatomical Pack with Tube	8.0 x 1.5 x 3.0 cm	Box of 10	10-120T
Anatomical Slimpack	8.0 x 1.0 x 3.0 cm	Box of 10	10-125
Anatomical Slimpack with Tube	8.0 x 1.0 x 3.0 cm	Box of 10	10-125T

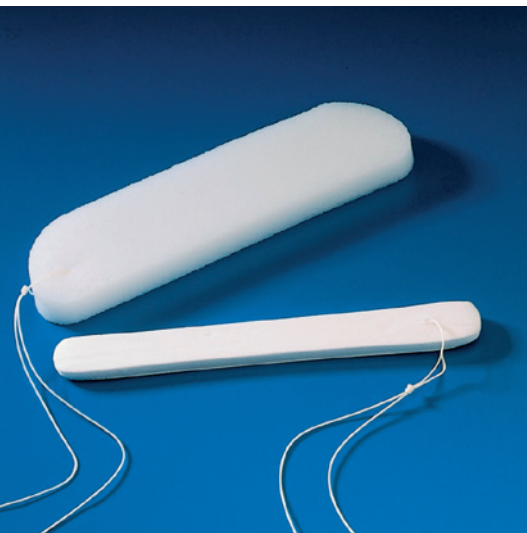

# PVA Epistaxis Nasal Pack

- Preoperative/emergency stenting tamponade
- Specifically designed to quickly contain and restrict blood loss
- Long (10cm) and Medium (8.0cm) Version to quickly pack entire nasal cavity
- Small (5.5cm) version to treat anterior epistaxis
- Layered Slimpack option allows thickness to be customised to individual requirements
- Suitable for delivering medication directly to the surgical area

Supplied sterile, single-use

Length: 5.5cm



Description	Dimensions	Pack Size	Code
Epistaxis Pack Short	5.5 x 1.5 x 2.5 cm	Box of 20	<b>10-150</b>

Length: 8.0cm



Description	Dimensions	Pack Size	Code
Epistaxis Pack Medium	8.0 x 1.5 x 2.5 cm	Box of 10	<b>10-145</b>
Epistaxis Slimpack Medium	8.0 x 1.0 x 2.5 cm	Box of 10	<b>10-146</b>

Length: 10.0cm



Description	Dimensions	Pack Size	Code
Epistaxis Pack Long	10.0 x 1.5 x 2.5 cm	Box of 10	<b>10-140</b>
Epistaxis Slimpack Long	10.0 x 1.0 x 2.5 cm	Box of 10	<b>10-141</b>
Epistaxis Slimpack Long Layered	10.0 x 1.0 x 2.5 cm	Box of 10	<b>10-141L</b>

# Series 5000™ PVA Nasal Pack

The NETCELL® Series 5000™ PVA Nasal Pack has been clinically proven to reduce trauma and improve patient comfort



Supplied sterile, single-use

Length: 6.0cm




Description	Dimensions	Pack Size	Code
Series 5000 Single Hydration Nasal Pack	6.0 x 1.5 x 2.0 cm	Box of 10	15-105

Length: 8.0cm




Description	Dimensions	Pack Size	Code
Series 5000 Single Hydration Nasal Pack	8.0 x 1.5 x 2.0 cm	Box of 10	15-110
Series 5000 Single Hydration Nasal Slimpack	8.0 x 1.0 x 2.0 cm	Box of 10	15-111

- Postoperative stenting tamponade
- Unique patent protected non-stick and soft outer surface prevents tissue ingrowth
- Specifically designed for use after septal, turbinate and rhinoplasty surgery
- Provides gentle, equal pressure within the nasal cavity when absorbing postoperative bleeding
- Remains hydrated for up to 7 days without rehydration
- Double compressed for easy insertion (50% slimmer than competitor products)
- Outer coating can be used in conjunction with water or oil based medication

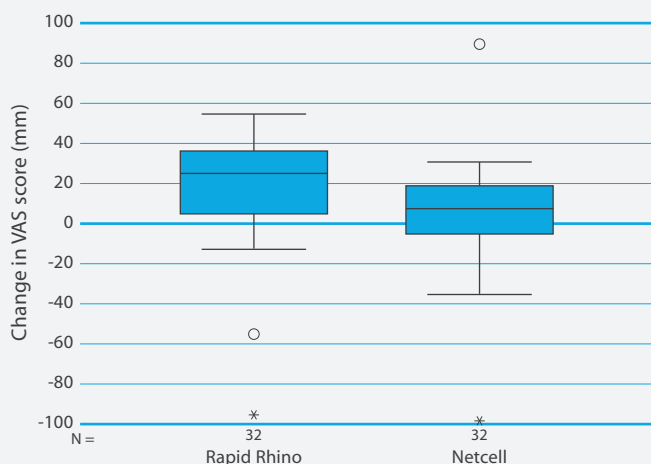

Box plots of change in VAS pain score for each pack type. The p value for the difference between these 2 data sets is 0.008.

\* Source International Rhinology 47, 4144, 2009.  
A randomised controlled trial comparing Rapid Rhino Mannheim and NETCELL Series 5000 packs following routine surgery

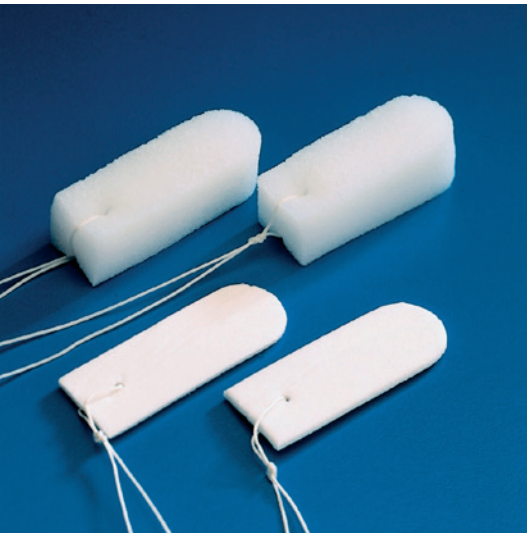

# PVA Sinus Pack

- Developed for postoperative packing following a F.E.S.S procedure
- Designed to prevent lateralisation of the middle turbinates postoperatively
- Can be used to pack the Osteomeatal complex to leave the nasal passage clear for easier breathing
- Suitable for delivering medication directly to the surgical area

Supplied sterile, single-use

Length: 2.5cm 

Description	Dimensions	Pack Size	Code
2.5 cm Sinus Pack Standard	2.5 x 1.2 x 2.0 cm	Pack of 2 / box of 10 packs	<b>20-200</b>

Length: 3.5cm 

Description	Dimensions	Pack Size	Code
3.5 cm Sinus Pack Standard	3.5 x 1.2 x 1.2 cm	Pack of 2 / box of 10	<b>20-210</b>
3.5 cm Sinus Pack Medium	3.5 x 0.9 x 1.2 cm	Pack of 2 / box of 10	<b>20-211</b>
3.5 cm Sinus Pack Small	3.5 x 0.6 x 1.2 cm	Pack of 2 / box of 10	<b>20-212</b>



# PVA Ear Wick

- Controlled expansion width-wise for optimal position and control
- Designed to facilitate instant delivery of antibiotics into the ear canal, to provide fast pain relief when treating Otitis Externa.
- Fenestrated option available to enable aeration



Supplied sterile, single-use

Length: 1.2cm 

Description	Dimensions	Pack Size	Code
12 mm Ear Wick	7 mm dia x 12 mm	Box of 10	30-301

Length: 1.5cm 

Description	Dimensions	Pack Size	Code
15 mm Ear wick	9 mm dia x 15 mm	Box of 50	30-300
15 mm Ear wick Fenestrated	9 mm dia x 15 mm	Box of 50	30-305

Length: 2.4cm 

Description	Dimensions	Pack Size	Code
24 mm Ear wick	9 mm dia x 24 mm	Box of 10	30-310
24 mm Ear wick Fenestrated	9 mm dia x 24 mm	Box of 10	30-315



## PVA Ear Packs

- Designed to prevent stenosis of the ear canal following surgery
- Controlled expansion width-wise for optimal position and control
- Can be used to support a graft following tympanoplasty procedure
- Available with string or without string

Supplied sterile, single-use

Length: 1.5cm 

Description	Dimensions	Pack Size	Code
Ear Pack (Strung)	9 mm dia x 15 mm	Box of 10	30-320

Length: 2.4cm 

Description	Dimensions	Pack Size	Code
Ear Pack (Unstrung)	12 mm dia x 24 mm	Box of 10	30-340



## PVA Microspheres

- Developed to control haemostasis during ENT microsurgery
- Also suitable for use during middle ear procedures to remove particles


Supplied sterile, single-use

Diameter: 3mm 

Description	Dimensions	Pack Size	Code
Microsphere Small	3 mm dia	Pack of 10 / box of 10	30-380

# PVA Instrument Wipes

- Specifically developed to effectively clean and extend the life of delicate microsurgical instruments
- Quickly removes protein, blood and particles without damage to the instrument
- Soft, flexible and strong material that is snag and fibre free

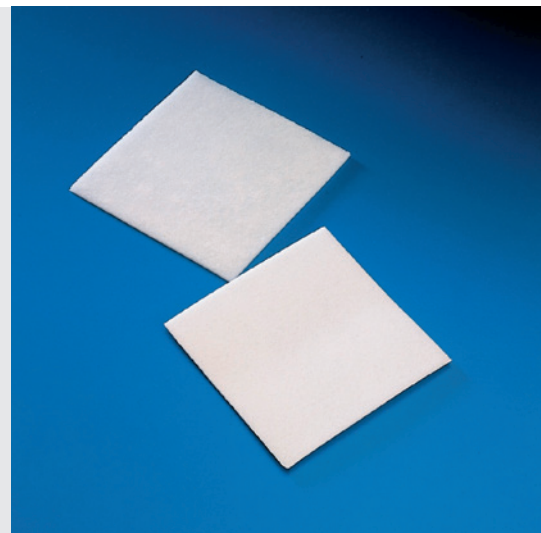

Supplied sterile, single-use

Description	Sterile / Non Sterile	Pack Size	Code
Instrument Wipe	Sterile	Pack of 1 / box of 20	35-900

# PVA Headband

- Beneficial for surgeons operating in demanding situations
- Provides an immediate cooling effect via evaporation once hydrated with water
- Reduces contamination to the operating site by controlling and limiting perspiration



Supplied sterile, single-use

Description	Sterile / Non Sterile	Pack Size	Code
Headband	Non-sterile	Pack of 1 / box of 20	35-950

**NETCELL®**  
ENT

**Network Medical Products, an Innovia Medical Company**  
Coronet House, Kearsley Road, Ripon, North Yorkshire, HG4 2SG, UK

**International** +44 (0) 1765 609555 | **US** 1.888.229.2875

[www.innoviamedical.com](http://www.innoviamedical.com) | [www.networkmedical.co.uk](http://www.networkmedical.co.uk)